

OWNER'S MANUAL

Model: F600B BARIATRIC FULL BODY PATIENT LIFT





Span-America Medical Systems, Inc. Greenville, SC 29606 (800) 888-6752

BEFORE USE:



Ensure the product has been assembled according to the instructions in this manual. All operators should read and understand the instructions for safe and proper operation of the patient lift.

P16790



Thank you for choosing Span

To better serve you, please record the following information for future use:

Supplier Name:	
Supplier Telephone:	
Product Serial Number:	
Date of Purchase:	

Date of Manufacture



The manufacture date can be found on the serial number. The YY on the serial number is the manufacture year, and the MM is the manufacture month.



This product has been tested by SGS, an accredited testing institute to comply with the requirements of:

ISO 10535, IEC 60601-1, ANSI/AAMI ES60601-1, CAN/CSA-C22.2 No.60601-1, IEC 60601-1-11, ANSI/AAMI HA60601-1-11, CAN/CSA-C22.2 No.60601-1-11, IEC 60601-1-2



Product Description

The Span-America F600B allow patients up to 600lbs to be lifted and transferred safely with minimal physical effort provided by the operator. That is due to it being a full electric lift. The F600B can lift the patient with ease from the floor to the bed. The 6 point spreader bar provide added comfort for the patient while moving from any surface With the addition of the optional battery pack, you will never have to be worried about your lift not being charged and ready for use.

i Intended Use

The Span F600B is intended to be used with patients who are completely immobile or totally dependent on the patient lift for transfers. The F600B is intended to be used for transfers of the patient in and out of bed, their wheelchair, to and from the commode, or any other type of surface. The F600B is intended to be operated by at least one caregiver, who will perform the lift and transfer. Transfers should always be done on a flat, dry surface and are not intended to be used in high humidity environments for extended amounts of time.

DEFINITIONS & SYMBOLS

In this manual the user refers to the patient or resident and may be used interchangeably at different times. Caregiver refers to the operator or person who is assisting with the transfer.

Symbols used in this manual and on the product and their meanings:

	-		
<u> </u>	Warning! Failure to heed this warning may result in damage to the product or serious injury to the operator and/or user.	×	Do Not Bleach. No bleach product may be used. The garment is not colorfast or structurally able to withstand any bleach.
i	Important instructions follow. Read and understand the instructions in the manual before using the product.	\otimes	Do Not Dryclean. Garment may not be commercially drycleaned.
†	Type BF applied part. Refer to the part of the medical device which come into physical contact with the patient in order for the device to carry out its intended	⋈	Do Not Iron. Item may not be smooted or finished with an iron.
	Class II symbol. To identify equipment is designed that a connection to electrical earth (ground) is not required.	\odot	Tumble Dry, Normal, Low Heat. A machine dryer may be regularly used at a maximum of Low Heat setting.
i	Note! Important information concerning the product and/or its correct and proper usage follows.	X	Do Not Dry. A machine dryer may not be used.
	Advise the operator to consult the accompanying documents before use.	180°F 82°C	Machine Wash, Hot. Initial water temperature should not exceed 82°C or 180°F.
11	Right side up	数	Do Not Wash
T	Fragile, handle with care	**	Keep dry
-25°C	Carton storage temperature limit of -25°C to 70°C	90%	Carton storage humidity limit of 0% to 90%
CE	European Conformity mark, The product conforms to European Medical Directive 93/42/EEC and meets applicable health, safety and environmental requirements.	**	Handle with care



Safety Instructions

A Span-America patient lift will provide years of service if it is properly maintained as any electric and/or mechanical equipment requires. Please pay careful attention to the following important information regarding the care, maintenance, and operation of the patient lift. Carefully read these instructions before assembling the lift, or attempting to lift a user with the device.

- Before attempting to lift any patient, please practice operating the lift and explain the lift procedures to the patient.
- Please provide special assistance under the caretakers' discretion.(i.e. patients who are comatose, spastic, agitated, or otherwise severely handicapped.)
- The patient lift should be used solely for transferring a user/patient from one utility (beds, chairs, toilets, etc.) to another. The patient lift should not be used for transporting or moving any patient from one location to another location.
- During lifting or lowering, whenever possible, always keep the base of the lift in the widest position.
- The base of the lift should be closed before moving the lift.
- Do not roll casters over any object while the user/patient is in the sling.
- · Do not lock casters during lifting.
- While being lifted in a sling, always keep the user/patient centered over the base and facing the caregiver operating the lifter.
- Never leave the user/patient unattended during lifting.

Service personnel may contact Span-America for any information, instructions, or certain parts for servicing purpose.



2




- **NO** modification of this equipment is allowed. It is unsafe to modify the equipment from its original state. This product is deemed safe to use with only the original parts as delivered. Unauthorized use of any modified equipment may result in serious injury or death.
- DO NOT replace any components of the lift without consulting with Span-America and must follow proper instruction from Span when replacing any components.
- **DO NOT** service any parts of the lift while in use with a patient.
- It is NOT safe to load the equipment more than the indicated maximum safe working load.
- IF any accessories with lower load capacity than the equipment itself is used,
 DO NOT load the equipment more than the lowest maximum safe working load.
- **REMOVE** the battery pack if the lift will not be used for some time, to prevent battery leakage. **PERIODICALLY** check the battery pack and the charger
- This equipment **DOES NOT** consist of any materials or ingredients to which patients are exposed, that may cause bodily harm or injury.
- Prolonged exposure or physical contact with any part of the equipment WILL NOT cause skin irritation, allergic reaction, etc.
- Using accessories, detachable parts, or materials not described in the instruction manual MAY RESULT IN SERIOUS INJURIES.
- Keep the equipment away from children and pets when in use. IT MAY CAUSE SERIOUS INJURY OR DEATH.
- Keep the equipment away from any pests or pest infested area. Diseases transferred from pests to the equipment MAY CAUSE SERIOUS ILLNESS to patients upon contact.
- If the lift is exposed to the sun for long periods of time, be CAUTIOUS as the surface of the lift may be hot.
- There are no risk if the lift is exposed to reasonably foreseeable environmental conditions
- **CAUTION!** The controller contains parts and assemblies that could be susceptible to **damage by electrostatic discharge (ESD)**. To reduce the likelihood of this, one should be grounded prior to using the controller.
- **DO NOT** use the lift on a inclined surface that is greater than 5°.

FEATURES & OVERVIEW

(i)

Prior to assembly, unpack and inspect all parts from the shipping carton. Contact you dealer immediately if any parts are damaged or missing.

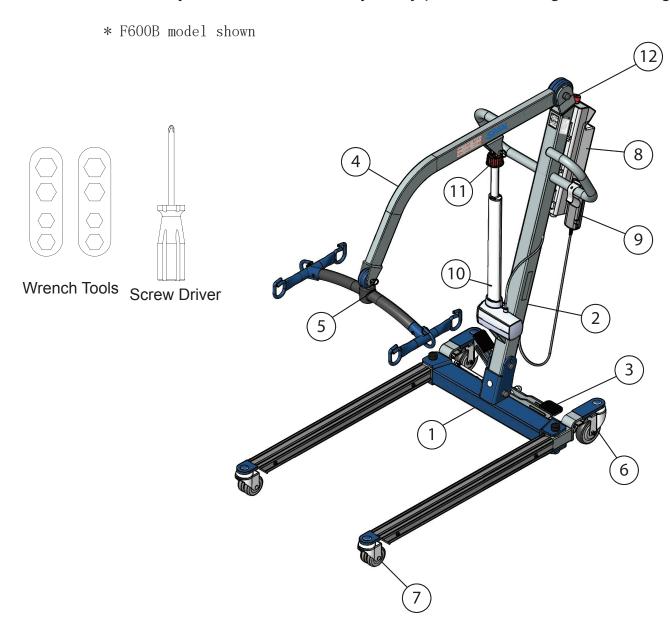


1	**Base	Combined	5	Spreader Bar	7.4 lb (3.4 kg)	9	Hand Control (P11620)	0.7 lb (0.3 kg)
2	Mast	Weight of 82 lb	6	Rear Caster (P16795)		10	Actuator (P16815)	10.2 lb (4.6 kg)
3	Leg Opening Pedal	(37.3 kg)	7	Front Caster (P16794)		11	Emergency Lowering	
4	Boom	14.6 lb (6.6 kg)	8	Control Unit (P16694)	11 lb (5 kg)	12	Emergency Stop Button	


^{**} Heaviest detachable part



Specifications of Control Unit Input Voltage 220-240VAC 50Hz / 110-120VAC 60Hz Output Voltage 24 VDC Battery Pack 24 VDC 5AH Lift Control Handset and Panel Protection Class IP54 (Protection against dust and splash of water) Sound Level Under Full Load 16 db Sound Level Under No Load 11 db Over-Duty Warning by LCD Battery Status by LCD

SPECIFICATIONS & OPTIONS

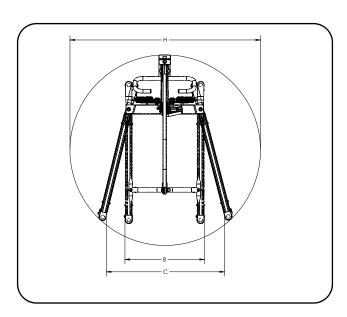

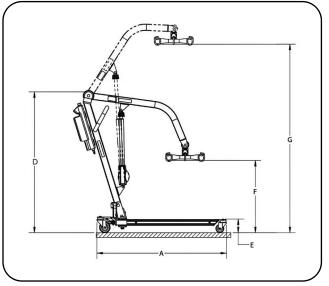



	Dimensions & Specific	ations
	Model	F600B
Safe Working Load		600 lb 273 kg
Α	Base Length	47.3" 1201.4 mm
В	Base Ext. Width - Closed	26.4" 670.6 mm
С	Base Int. Width - Opened	39.2" 995.7 mm
D	Mast Height	51.3" 1303.0 mm
Е	Front Caster Height	4.7" 119.4 mm
F	Lifting Height - Lowest	26.3" 668.0 mm
G	Lifting Height - Highest	68.7" 1745.0 mm
Н	Turning Diameter	63.4" 1610.4 mm

Sta	ndard Features
1	6-Point Spreader Bar
2	4" Rear Caster with Brake
3	3" Dual Front Caster
4	Power Control Unit
5	Reinforced Foot Pedal
Op	tional Features
1	Charging Cradle
2	Digital Scale (P11634)
3	Additional Battery pack

	F600B
Lift Weight Unloaded (lb)	123
Lift Weight Unloaded (kg)	56

^{*} Span-America is committed to continuous improvements of our products therefore the specification, dimensions, and features listed above are for guidance only and are subject to change without prior notice.

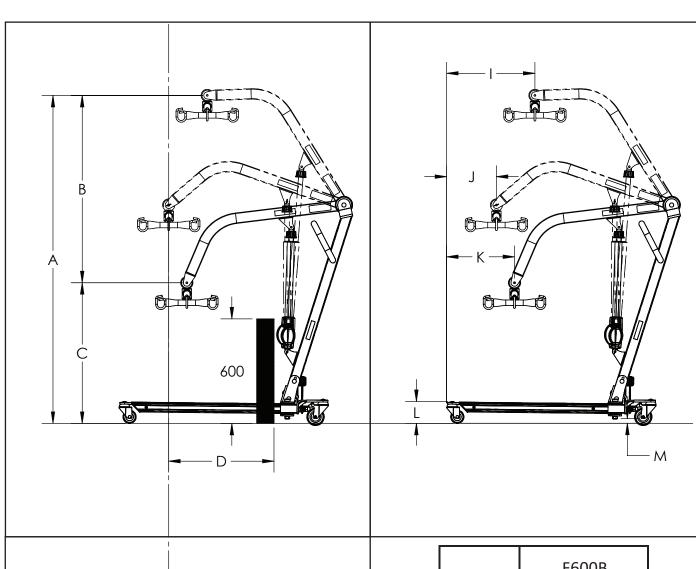

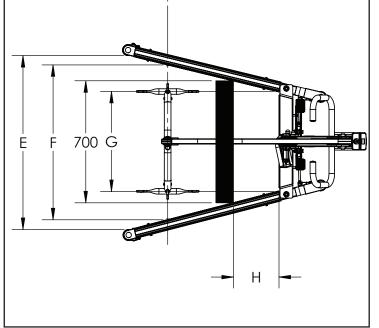


	F600B	
	mm (in)	
А	1882.4 (74.1)	
В	1074.9 (42.3)	
С	807.5 (31.8)	
D	604.8 (23.8)	
Е	1061.5 (41.8)	
F	889.0 (35.0)	
G	574.6 (22.6)	
Н	262.0 (10.3)	
I	506.5 (19.9)	
J	285.5 (11.2)	
K	389.7 (15.3)	
L	126.0 (5.0)	
М	26.2 (1.0)	

6


/

F600B Installation


Step 1: Remove folded lift from box.

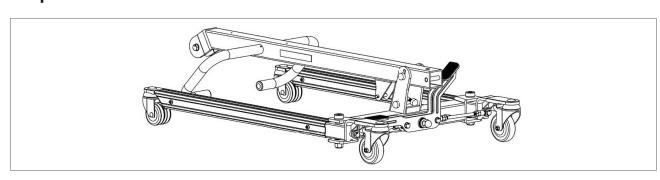


Step 2: Remove the bottom plastic caps, then the bottom bolt, nut and washers on the side of the base bracket. Use the wrench tool to loosen the bolt if needed.

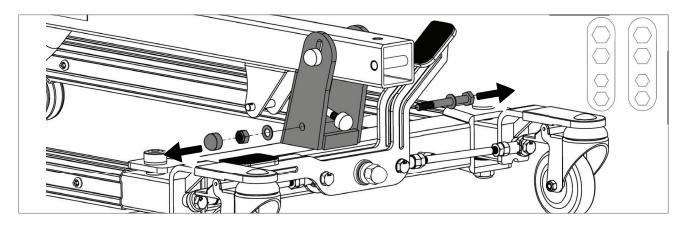


Step 3: Remove the plastic caps, then the bolt, nut and washers at the front of the base bracket. Use the wrench tools if needed.

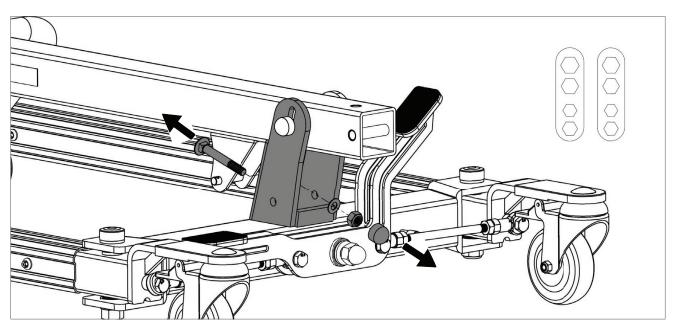


Step 4: To stand the mast upright, lift the mast upward 1" to clear the metal plate at the bottom of the mast.

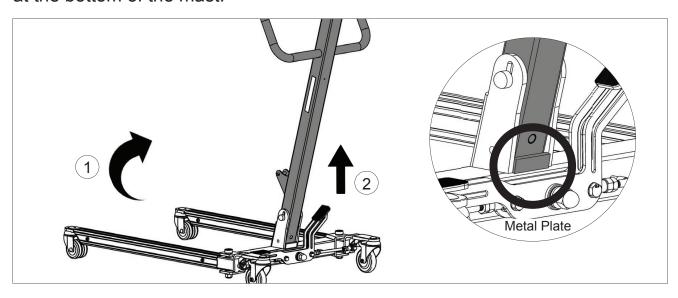


Step 5: Reinsert and tighten the front and side bolts, nuts, washers to secure the mast using the wrench tools. Reinsert the plastic caps after the screws are tightened.

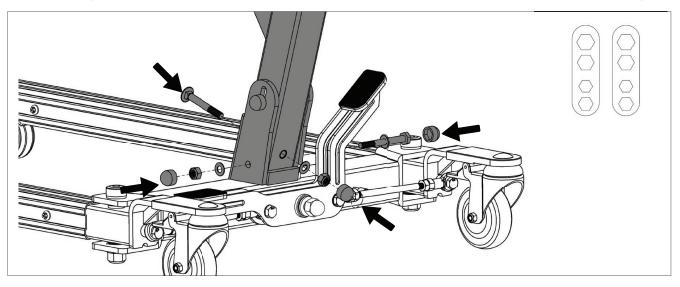


Step 6: Remove the plastic cap, and then bolt, nut and washer at the top of the mast. Use wrench tools if needed.

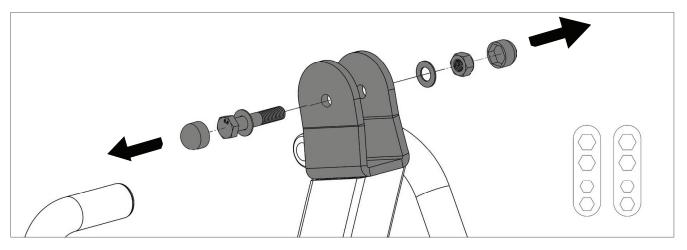


Step 7: Install the round plastic cover and plastic washers to the boom.

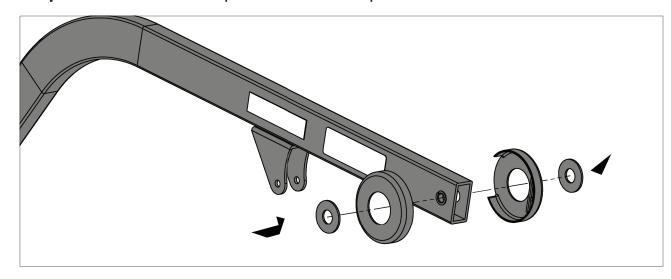


Step 8: Install the boom to the mast and align the holes. Then secure it by reinstalling the bolt, washer and nut. Once the screw is tightened using the wrench tools, reinsert the plastic cap.

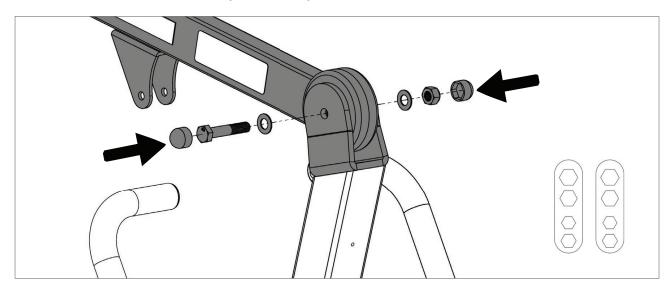


Step 9: Secure the lower end of the actuator to the bracket on the mast using a pin and retaining ring as follows.

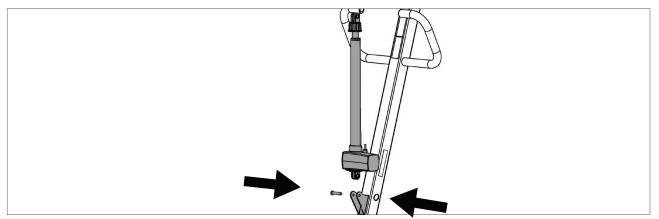


Step 10: Secure the upper end of the actuator to the bracket on the boom using a pin and retaining ring as follows.

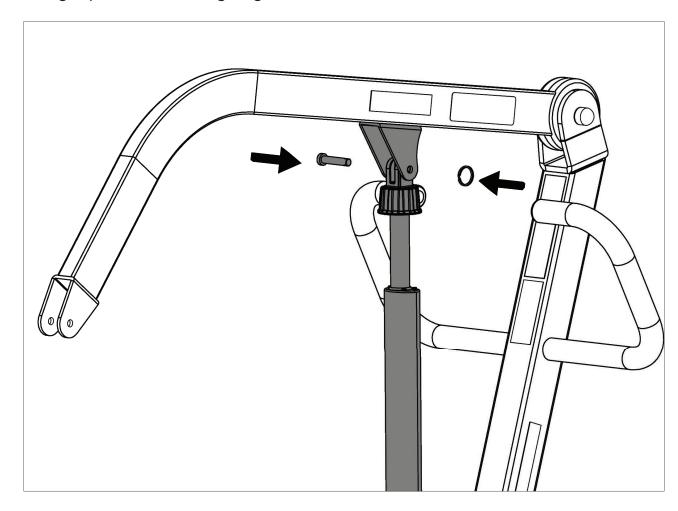

Step 11: Secure the spreader bar to the boom using a pin and retaining ring as follows.

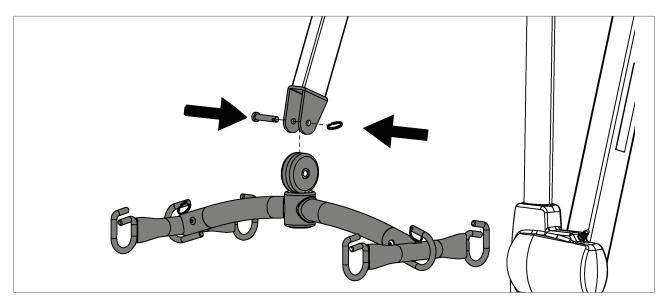

Power Control Unit Installation

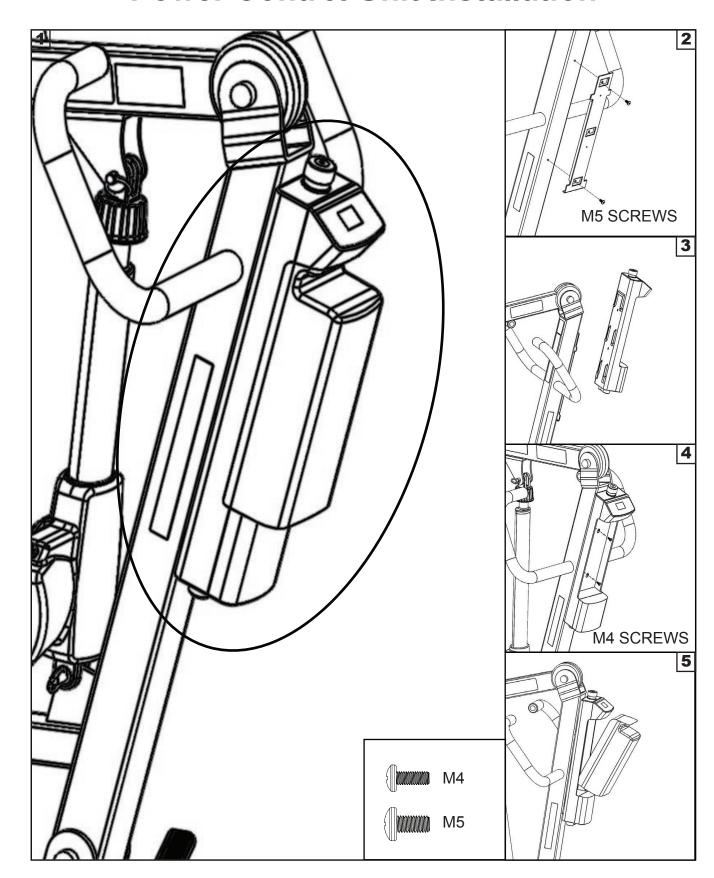

Operating Instructions





Double check all assemblies for tightness and read operating instructions carefully prior to use.

For optimum performance the lift should be transported and stored in following condition range:

- -25°C to +5°C (-13°F to 41°F), and
- +5°C to +35°C (41°F to 95°F) at a non-condensing relative humidity 0% to 90%
- +35°C to +70°C at a water vapor pressure up to 5 kPa

For optimum performance the lift should be operated in following condition range:

- at a temperature range of +5°C to +40°C (+41°F to +104°F)
- at a non-condensing relative humidity range of 15 % to 90 %, but not requiring a water vapor partial pressure greater than 5 kPa
- 70 kPa to 106 kPa atmospheric pressure

The lift will not have any potential electromagnetic or other interference with other devices.




PREPARATION BEFORE LIFTING

- Twist the RED EMERGENCY STOP BUTTON clockwise.
- Verify the Lifting and Lowering Functions by pressing either the UP or DOWN button on the hand control once.
- Position the base of the lift around or under the object.
- Widen the base and disengage the caster brakes.



Do not attempt to transfer a patient without prior approval of the person's nurse. Also, do not transfer without having studied the instructions and performed several practices in operating the product. Together (with the patients doctor, nurse, or medical attendant) select a sling that is both practical and comfortable. The sling selected should be one that serves the needs of the patient, while providing the patient with optimal safety. Never interfere with the lift, unless instructed by the attendant. Have a doctor, nurse, or medical attendant present during the first few times the lift is used to transfer a new user.

Power Control Unit Overview

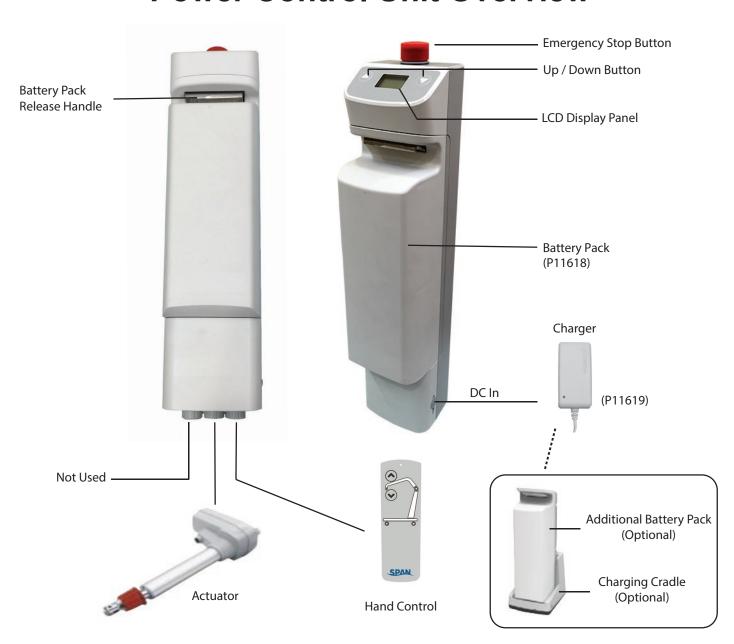

- 1. Connect Actuator as shown above.
- 2. Connect Hand Control as shown above.

<u>^</u>

WARNING!

- Contains a long charger cord that may cause strangulation if used improperly
- Keep long charger cord out of walkway where it poses a tripping hazard
- Have the charger port on the controller facing away from the wall for easy disconnection






LCD Display Panel Signs

When the emergency stop button is released, the LCD Display Panel will show one of the four signs below.

- The sign will display for 5 seconds.
- Then the lift will go into standby mode and the sign disappears.
- The lift is ready for use.



Charging: DC Charger or Charging Cradle

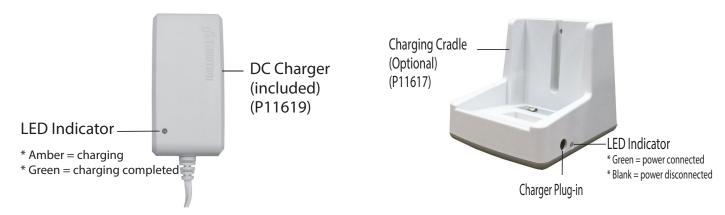

LCD Display Panel Signs - Charging

During charging, the LCD Display Panel will show one of the four signs below.

- The DC Charger must be connected.
- The emergency stop button must be released.
- The battery will not be charged if the emergency button is pressed in.




Operating Lift: Using Hand Control




Button	Press Button	Release Button
•	▼ ∴	

Operating Lift: Using LCD Display Panel

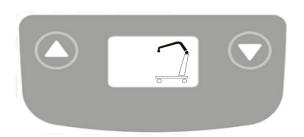

Button	Press Button	Release Button
	♣	


Warning! - Battery Low and Charging is Needed

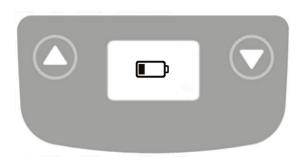

If the battery needs to be charged, the LCD Display Panel will show a blinking low battery sign shown on the left either when the emergency stop button is released or when a button on the LCD Display Panel is pressed.

- Lift is still functioning.
- Buzzer will beep every 2 seconds.
- Battery needs to be charged immediately after use.


Warning! - Overloaded

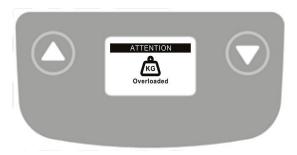

During lifting if the lift is overloaded, the LCD Display Panel will show the sign on the left and lifting will be deactivated immediately.

- The buzzer will beep 3 times.
- Lowering the lift by pressing the DOWN button is still allowed.
- The emergency stop button must be pressed in and then release the button before further lifting.

Warning! - Charge Battery Immediately

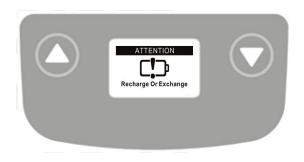

If the battery is critically low, the LCD Display Panel will show the sign on the left either when the emergency stop button is released or when a button on the LCD Display Panel is pressed.

- Lift operation is deactivated.
- Battery needs to be charged immediately before further use.

16

Available slings compatible with Patient Lift

Sling Type	Material	Sizes	
Universal with head support (Span-America)	Mesh with padding	S (M920276) XS (M920930) M (M920277) XL (M920931) L (M920278) XXL (M920932)	
Universal without head support (Span-America)	Mesh with padding	S (M920667) XS (M920936) M (M920668) XL (M920937) L (M920669) XXL (M920938)	
Hygienic Sling (Span-America)	Polyester	S (M920279) XS (M920939) M (M920280) XL (M920940) L (M920281) XXL (M920941)	
Disposable universal with head support (Span-America)	Nonwoven, Disposable	S (M920942) XL (M920945) M (M920943) XXL (M920946) L (M920944)	
Sani	Polyester	S, M, L, XXL	
Invacare® Replacement	Polyester	M, L, XL	
mvacare* kepiacement	Mesh	M, L, XXL	
Invacare® Replacement with Head Support	Mesh	M, L, XL, XXL	
Hoyer® Replacement	Polyester	M, L, XL	
Hoyer® Replacement with Head Support	Polyester	M, L, XL	
Liko®/Guldmann® Replacement	Polyester with PU Coating	M, L, XL	
Liko®/Guldmann® Replacement with Head Support	Polyester with PU Coating	M, L, XL	
Liko®/Guldmann® Disposable Replacement	Nonwoven, Disposable	M, L, XL	
Liko®/Guldmann® Disposable Replacement with Head Support	Nonwoven, Disposable	M, L, XL	
Dodiatria with Hood Cupport	Mesh	Ctandard	
Pediatric with Head Support	Polyester	Standard	



If the patient is pregnant or has skin conditions, consult with a doctor before use.




ALWAYS check that the sling weight capacity is suitable for the patient



All slings are type BF applied parts

Sling Hook-Up Guide

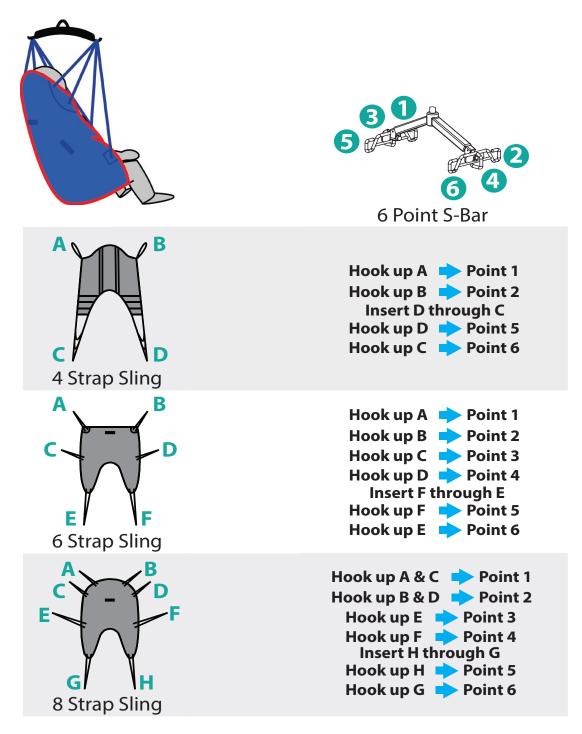



BEFORE LIFTING THE PATIENT, PULL ON **EACH** SLING LOOPS TO ENSURE THE LOOP IS HOOKED ONTO THE SPREADER BAR PROPERLY

For further product information and support, please contact your dealer or visit our website at www.spanamerica.com

(

Transfer From Bed to Wheelchair

- 1. User/patient should be in the center of the bed.
- 2. Position user onto his/her side by rolling user towards you.
- 3. Roll the sling in half. The handle on the back section sling should face outward when the sling is fitted.
- 4. Position the sling under the user/patient so the commode aperture aligns with the base of the spine and top of the sling close to the neck.
- 5. Roll user onto the opposite side and position him/her on the flat section of the sling. This will allow you to unroll the remainder of the sling from the other side of the user.
- 6. Once the sling has been positioned in the center of the user/patient, feed the leg sections under the thighs and draw them up between the thighs.
- 7. Raise the head of the bed if this function is available.
- 8. Move the lift slowly towards the user and position the spreader bar over the user's/patient's chest. Attach clips of the sling to the hooks on Spreader Bar as shown in the Sling Hook-Up Guide.
- 9. Lift the user above the bed by using the hand control.
- 10. Pull lift away from bed. Position user over the wheelchair or chair then lower the patient onto the surface.
- 11. Reverse the above procedures when returning the user/patient to bed.

Transfer From Wheelchair to Bed

- 1. Help the user lean forward slightly, then slide the sling down between the chair and the user's back.
- 2. Position the commode aperture where the buttocks meets the seat.
- 3. Position the sling equally around both sides of the body.
- 4. Draw the leg sections to the front along the length of the user's thigh.
- 5. Check the sling's central positioning by comparing the lengths of the leg sections when they are drawn forward.
- 6. Reposition the sling if the leg sections are not equal in length.
- 7. Feed the leg sections under the thighs.
- 8. From between the legs, gently pull the leg section up the inner thigh. Feed as much material as possible under and between the thighs. Ensure the leg sections are positioned midway under the thighs to provide good support and greater comfort.
- 9. Move the lift slowly towards the user and position the spreader bar over the user's chest. Attach clips of the sling to the hooks on Spreader Bar as shown in the Sling Hook-Up Guide.
- 10. Lift patient above the wheelchair by using the hand control.
- 11. Pull lift away from wheelchair. Position patient over bed and lower patient onto bed
- 12. Reverse the above procedures when returning the user/patient to a wheelchair.

Maintenance & Inspection



- The operator of the lift shall inspect the Span lift before each use. Check all bolts for tightness. Make sure the base can be easily widened, and that all lift parts are in place. Make sure that casters can be turned freely, and that caster brakes can be engaged. Make certain all necessary items (i.e. slings, and wheelchairs, etc.) are ready for use.
- At least once a month, the lift should be thoroughly inspected by a person qualified to recognize any signs of wear and tear, and looseness of bolts or parts. Replace any worn parts immediately.
- To lubricate, put silicone lubricating grease on the following points when the Span lift is placed into service and every month thereafter: Top of Mast, Spreader Hinge, and Caster Axles.
- With proper use and care, the expected lifetime of the lift is 10 years or 20,000 cycles.
- The expected lifetime of the electrical components are 3 years.

	First Received	Monthly	Every 3 Months
Boom & Spreader Bar			
Check connections between boom and		_	
spreader bar for improper connection,		•	
Check connections between boom and mast			
for improper connection, looseness or wear.	•		
Check the boom for bending and deflection.		•	
Mast			
Check the mast for bending or deflection.	•	•	
Base & Foot Pedal			
Check bolts and nuts for looseness.	•	•	
Check casters and axle bolts for tightness.	•	•	
Check rubber parts on the casters for			
deflection.	•		
Apply grease to caster ball bearings if needed.			•
Check welding joints for rust and crack.		•	•
Actuator & Control Box			
Make sure the control box is firmly affixed to		_	
the mast.		•	
Make sure the actuator is secured to the			
boom and mast with pins and key rings.	•		
Make sure actuator plug into control box is			
not loose.			
Sling & Sling Hardware			
Check sling for wear.	•	•	
Check sling hardware every time before use.	•	•	
Lift Cleaning	As Needed.	As Needed.	As Needed.



Cleaning and Disinfecting

Cleaning of the lift

- Use pH-neutral detergents only.
- Remove the battery from the controller
- Clean lift surfaces, control box and hand control with a damp cloth using pH-neutral detergent.
- Remove debris and hair from the casters and make sure the casters can spin and swivel freely.



Solvents or strong liquids or abrasive cleaning products will damage the surfaces of the lift.

Cleaning of the sling

- Check with sling label for washing condition.
- · As for minor cleaning between patient uses, we recommend wipes with water-soluble soap or isopropyl alcohol.
- For slings with rigid supports, remove the plastic inserts before washing. (Look for an opening near the plastic inserts for removal)
- Machine wash warm or cold.
- Do not wash with other colors.
- Air dry or very low dryer heat.



- · DO NOT USE BLEACH.
- DO NOT clean or wash disposable slings.REPLACE a sling when it show signs of deterioration.

Deterioration of sling

· Contact Span-America for "Guideline for Identifying Deteriorated Slings" if needed.

Disinfecting of the lift and sling

- Inspect slings prior to each use for contamination from previous use.
- Use isopropyl alcohol for any disinfection purpose.




- Wash a reusable sling for serious contamination
- Replace a disposable sling for any contamination
- Reusing a contaminated sling can cause serious bacterial infection

Contact your dealer immediately if standard troubleshooting techniques do not correct the failure. Do not attempt to lift until all failure and safety issues have been resolved.

In case of lift failure, please follow the procedures below to safely lower the user.

The Emergency Lowering Device is located at the top of the actuator shaft. It is intended for use if the actuator or electronics fail to operate while the patient is suspended in midair.

The device consists of a plastic collar ring that should be first pulled down and turned clock-wise continually until the patient has been lowered.

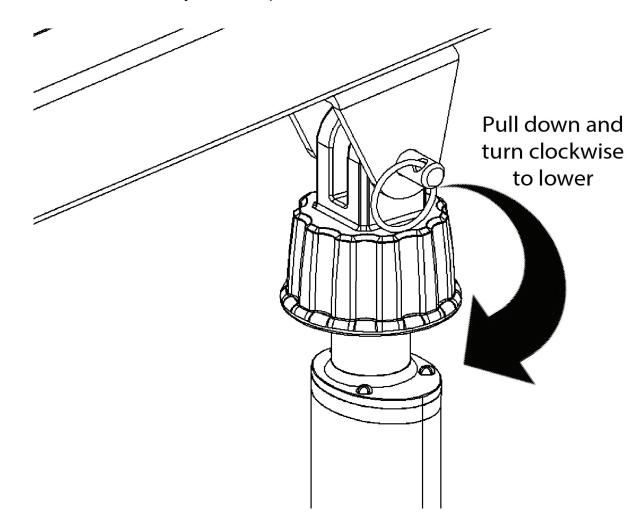

NOTE:

THE EMERGENCY LOWERING DEVICE IS INTENDED FOR USE DURING LIFT FAILURE. THIS DEVICE WILL ALLOW LOWERING OF PATIENTS ONLY. PLEASE CONTACT YOUR DEALER IMMEDIATELY IN CASE OF FAILURE.

Troubleshooting Guide

The following list of encountered problems and solutions will assist you in determining what may be causing your patient lift not to function as designed. If you have a problem occurring which is not listed below please contact your dealer or technical support for help. Do not attempt to repair or replace components or parts on your lift as this may void your warranty or cause further problems that may result in patient injury. Stop using your lift immediately if it is not functioning correctly or any warning beeps are heard.

I need to re-charge my battery pack often or it fails to hold a charge when I charge it.

- Replace the battery pack as it is at the end of its life cycle
- Battery pack should be changed every 18-24 months depending on usage

The actuator will either go up or down but not both

- Replace your hand control as it is at the end of its life cycle
- Your control box may be malfunctioning

There is a grinding sound inside the actuator when lifting

Replace your actuator as the internal gears are stripped

The actuator stops and starts while lifting or lowering

Replace the batteries or battery pack as they are at the end of their life cycle

My lift will not operate even though it shows the batteries are charged

- Unplug the AC power cord from the control box as the lift will not work when the control box/ charger is connected to AC Power.
- Check to make sure the hand control is properly inserted all the way into the control box port.
- Check to make sure the actuator plug is properly connected and inserted all the way into the control box
- Your hand control needs to be replaced as it is at the end of its life cycle

If there is a problem you could not solve, please contact your dealer for assistance.



WARRANTY POLICY

Span-America ("Span") offers a limited warranty on all patient lifts, slings and accessories to be free of defects in workmanship and product performance. This warranty extends only to the original purchaser and is non-transferrable. All warranty claims must be submitted by the authorized dealer or distributor who originally sold the product with proper proof of sale and serial number where applicable. Prior Return Authorization (RA) from Span is required for all warranty replacements. Span reserves the right to repair or replace only defective parts or accessories in lieu of a complete new patient lift. The repaired or replacement part shall be warranted for a period equal to the remainder of the warranty period of the defective part.

Span does not provide advance replacements for warranty claims. For situations where an end user requires a replacement in advance and before a warranty item can be returned, dealer or distributor must purchase the replacement at full cost. All freight charges for the replacement are the sole responsibility of the dealer or distributor. When end user returns warranty item and item proves to be defective, Span will issue full credit for the replacement less freight charges.

Product returned without RA number clearly marked on the package or product returned later than 30 days after authorization will be refused and returned at sender's expense. Final disposition of warranty claims will be determined at the sole discretion of Span. Warranty claims will be denied for any of the following: product abuse or misuse, accidental or malicious damage, improper installation, product used with parts, components or accessories with quality or specifications incompatible with product, adulterated product, user neglect, failure to maintain and service product as specified in the owner's manual or care tag, serial number removed or defaced, or normal wear and tear.

Span warrants the following products and components beginning from the purchase invoice date for the specified time period:

WARRANTY PERIOD	PRODUCT OR COMPONENT
Three (3) Years	Patient lift frame or spreader bar
Two (2) Years	Actuator, control box, pendant, charger, weigh scale, or casters excluding normal tread wear
One (1) Year	Battery pack
Six (6) Months	Reusable fabric slings

Single Patient Specific slings are designed for limited use with one patient and may not be laundered. Span will replace any disposable sling found to have a manufacturing defect. Normal wear and tear will not be covered under warranty.

This warranty expressly excludes wearable components including but not limited to foam parts and caster tread. Furthermore, this warranty is void and null for product that has not been purchased or paid for in full.

RETURN GOODS POLICY

Patient lifts may not be returned unless the wrong lift is shipped in error by Span or the lift is heavily damaged or defective out of the box. For all other items, purchaser may request a RA for purchased goods within thirty (30) days of purchase invoice date. All returns must be received by Span no later than thirty (30) days after authorization or the RA will be voided. Return package must be clearly marked with the RA number or the package may be refused and returned at sender's expense. Patient lifts are subject to a minimum restocking fee of 25% or more. Please note that patient lifts being returned must be in the original carton with all parts, components and packing materials included. Failure to comply with this requirement will incur higher restocking fees or a rejection of the return. Slings, parts and accessories may only be returned if they have not been used. There are no exceptions to this provision.

All freight charges are the sole responsibility of purchaser when any of the following occurs: (i) ordering error where an incorrect item is shipped in accordance with purchase order, (ii) an order is cancelled while in transit or (iii) delivery is refused by customer.

Span reserves the right to issue credit amounts based on strict adherence to this policy.









EC REP

Advena, Ltd. Tower Business Centre, 2nd Flr, Tower Street, Swatar, BKR 4013, Malta

Span-America Medical Systems, Inc. Greenville, SC 29606 (800) 888-6752

P16790

Australian Sponsor: KD&A Pty Ltd. 286 Flinders Street Adelaide, South Australia, 5000 (ACN 114 941 338)

ECO # 2329

www.spanamerica.com | © 2020 Span-America

